


# ClearVial™ Ready To Use – packed in PET trays

Substitution for glass

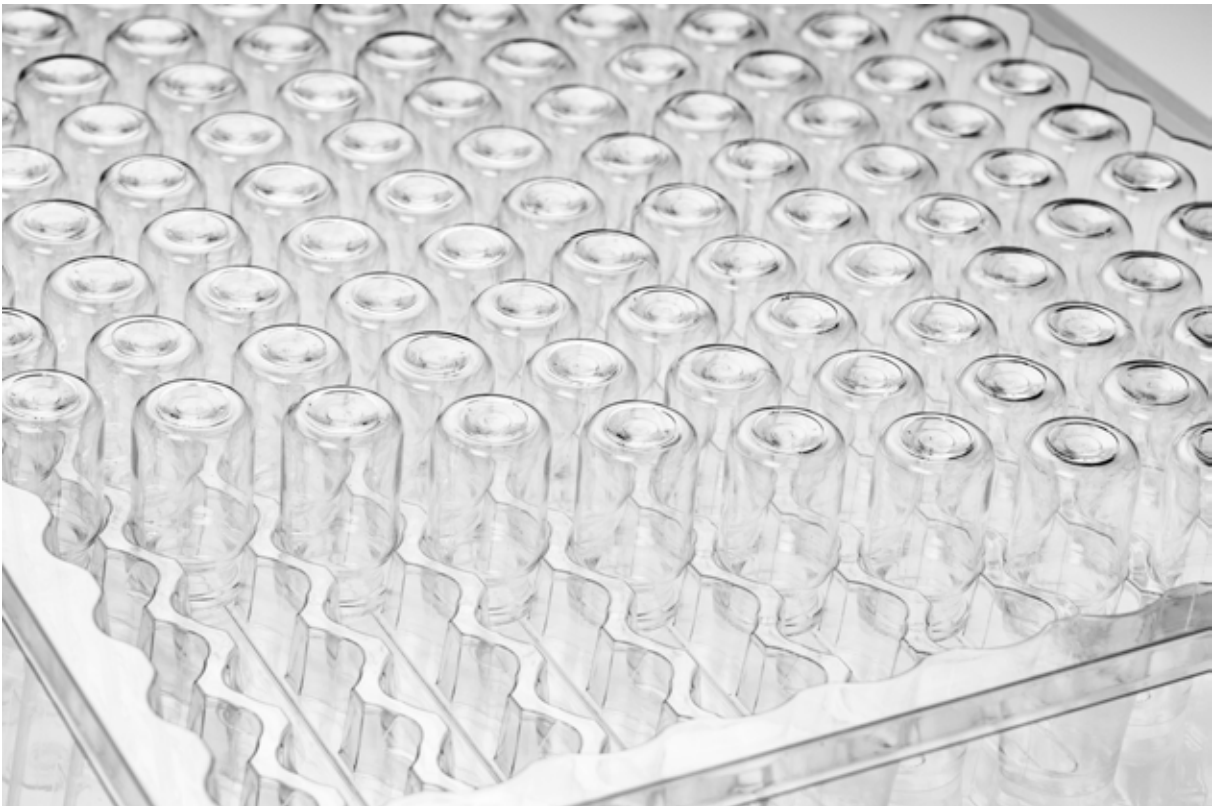

ClearVial™ Ready To Use – Packed in PET trays

## **MedicoPack introduces new packaging method for COC vials**

For more than 10 years, MedicoPack has successfully supplied COC vials packed in shrink and vacuum sealed pouches to the global pharmaceutical industry. But in line with an increasing demand for tray solutions, we have chosen to develop a new packaging method, while also taking customers' requirements into account.

That's why MedicoPack is proud to introduce our brand new standard PET tray solution.

# Clear Vial™

**ClearVial™ made of COC is a reliable and cost-effective alternative to glass vials.**

Transparent, lightweight and with good impact strength, ClearVial™ is a breakthrough in primary packaging solutions for the pharmaceutical and biotech industries. Using the amorphous polymer, COC (Cyclic Olefin Copolymer), transparency is no longer an issue with plastic vials. ClearVial™ is delivered sterile and can be used directly in the filling line with no need for prewash, drying and sterilization. This saves water and energy. ClearVial™ comes in a standard range of 5 mL, 10 mL, 20 mL, 50 mL and 100 mL following ISO 1862 standard.



## **Benefits of ClearVial™**

- Packed in clean room ISO class 7 / EU-GMP class C.
- Supplied Ready To Use - no wash needed before filling.
- Is 100% transparent, enabling easy inspection.
- Offers good barrier properties.
- High impact strength - provides a safer work environment.

# Vials packed in PET trays



## Benefits of packing in trays

- No vial to vial contact minimizes the risk of scratches and transport damages.
- Outer dimensions identical for all tray sizes -> easy handling at your filling line.
- Vials ready for pre-sterilization by irradiation or EtO.
- PET trays can be recycled as rPET for beverage packaging.

## Easy to peel

The lid is made of Tyvek®, characterized by

- Good seal strength.
- Easy opening.
- Excellent peelability.



## Packing configuration

Product number	1875T	1876T	1877T	1879T	1880T
Volume (mL)	5	10	20	50	100
Vials per tray (pcs)	160	149	98	50	30
Vials per carton (pcs)	2,560	2,384	1,372	500	240
Quantity per EUR pallet (pcs)	30,720	28,608	16,464	6,000	2,880

The data set forth herein has been carefully compiled by MedicoPack A/S. However, there is no warranty of any kind either express or implied, applicable to its use and the user assumes all risk and liability in connection therewith.

